

Usage

01. Before use, please check if the product can be started. If it cannot be started, please first charge the product with the included charging cable before use.
02. Attach the massager to the knee or waist pads,
03. Press and hold the power button for 3 seconds, and the buzzer will beep once. The default temperature for power on is 45%, and the default massage intensity is set to level 01.
04. When the product is turned on, the default temperature is 45 °C. Press the "cycle button" once to increase the temperature, and the temperature will cycle in multiple levels. The default massage intensity is Press the "cycle button" in gear 01 to increase the massage intensity, and the massage intensity will cycle through multiple levels.
05. In the startup state, long press the power button, the buzzer beeps once, the product shuts down, and all indicator lights are turned off.

Product Introduction


Product Name: Massage Equipment

Rated voltage: 5V

Rated power: 9W

Battery capacity:

Note: This product is rechargeable and compatible with TYPE-C interface charging

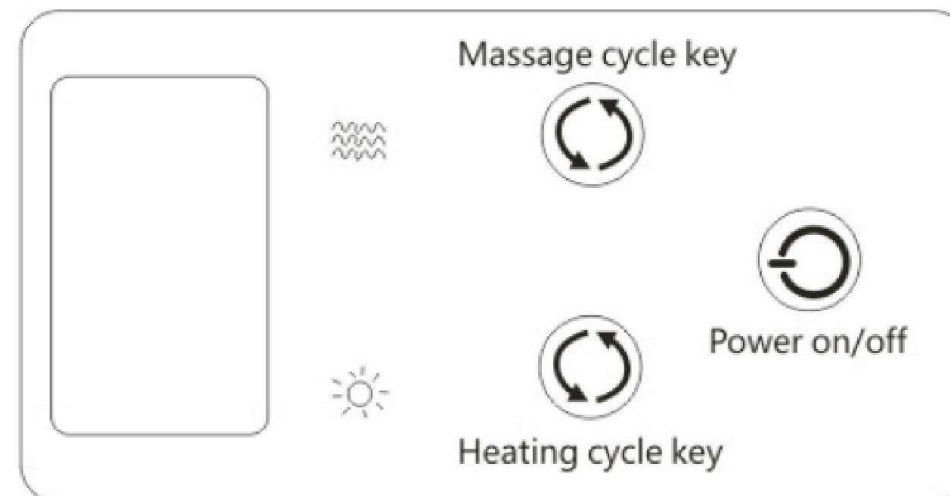

work environment

	work environment	Storage and transportation
Temperature range, humidity range, atmospheric pressure	0°C~40°C < 80%RH 70KPA-106KPA	-0°C~50°C < 80%RH 70KPA-106KPA